


Microbiology/Histology

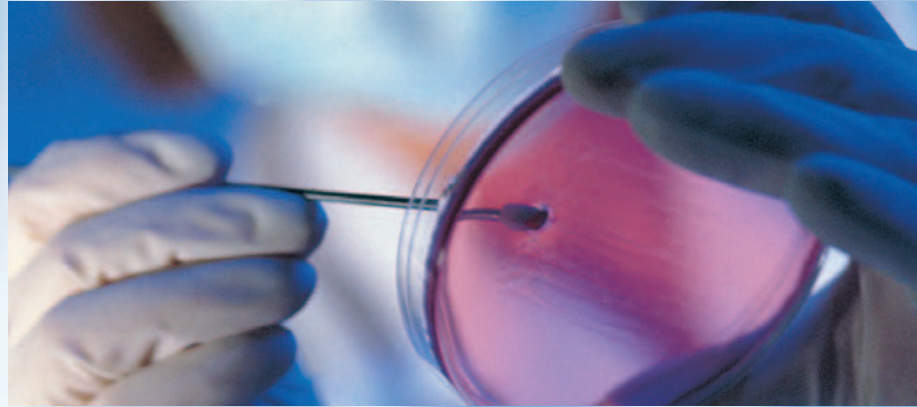

Loops

Needles

Streakers

Nichrome

Platinum

Clearing Agents

Formalin Neutralizer

Odor Control

Cytology Fixative

Freezing Media

www.deconlabs.com





CitriSolv® Solvent and Clearing Agent

CitriSolv is a d-Limonene-based solvent and clearing agent that can be used as a safe alternative to xylene and ethyl acetate. It is soluble with absolute alcohol and will dissolve paraffin, glues, adhesives, and mounting media. CitriSolv provides for better tissue sections – it will not harden or shrink tissue. It dries slowly with no residue and enhances some stains. Use in manual or automatic

tissue processor staining procedures in histology/cytology applications that require xylene.

Item	Standard Size	Alternate Size	Decon Catalog #
CitriSolv	1 gallon (3.8L)	4 x 1 gallon (3.8L)	1601
CitriSolv	5 gallon (20L)	N/A	1603



CitriSolv® Hybrid Solvent and Clearing Agent

CitriSolv Hybrid is a mixture of d-limonene and hydrocarbon material allowing for the reduction of orange aroma and an increase in flash point. Not classified as flammable.

Item	Standard Size	Alternate Size	Decon Catalog #
CitriSolv Hybrid	1 gallon (3.8L)	4 x 1 gallon (3.8L)	1601H



SeraSpin™ Crystals

Decon's SeraSpin Crystals are pellets that are used when centrifuging whole blood to help trap fibrin particles and clots in serum in order to maximize serum yield. Using these crystals eliminates the need to centrifuge a second time.

Item	Standard Size	Alternate Size	Decon Catalog #
SeraSpin	Net WT. 16 oz.	6 x 16 oz.	2116



Odorgone® & The Pill™ Autoclave Deodorizer

Heat and moisture from the autoclave release the contents of the capsule to over-

come unpleasant smells generated from heating microbiological waste. Simply place one capsule in each load to absorb offensive odors leaving the laboratory smelling fresh and clean. Odorgone comes in two pleasant scents; Lemon and Apple. The PILL is fresh scented.

Item	Standard Size	Alternate Size	Decon Catalog #
OdorGone–Apple	25 Capsules	8 x 25 Capsules	APP25
OdorGone–Apple	100 Capsules	6 x 100 Capsules	APP100
OdorGone–Lemon	25 Capsules	8 x 25 Capsules	LEM25
OdorGone–Lemon	100 Capsules	6 x 100 Capsules	LEM100
The PILL	100 Capsules	20 x 100 Capsules	5100
The PILL	500 Capsules	6 x 500 Capsules	5500



SlideRite™ Cytology Fixative

SlideRite is a highly concentrated, water-soluble spray that forms a lasting protective coating over all types of cytological smears. Many cytological specimens are fragile and degenerate quickly, therefore, an immediate fixation of the material to be studied is essential.

Item	Standard Size	Alternate Size	Decon Catalog #
SlideRite	1.5 oz.	12 x 1.5 oz.	1801



Flash Freeze™

Flash Freeze is moisture-free propelled cold gas used to rapidly cool or freeze samples in clinical or scientific laboratories.

Item	Standard Size	Alternate Size	Decon Catalog #
Flash Freeze	11 oz.	6 x 11 oz.	3511



Inoculating Loops, Streakers and Needles

Made of high quality, twisted nickel chromium wire, LeLoops and LeStreaks are re-usable making them more economical

than disposable plastic loops. Loops slide smoothly across the agar surface without tearing or gouging. The nichrome wire cools rapidly after flaming, saving time. Loops are available calibrated or non-calibrated. All calibrated loops are subjected to Evan's Blue dye test.

LeLoops come without holders, which are available separately. Holders are made of nickel-plated brass with well insulated handles.

LeStreak and LeStab are permanently attached to lightweight, insulated aluminum holders perfectly balanced for easy manipulation.

LeLoops®

Item	Loop ID	Std. Size	Alt. Size	Decon Cat #
Calibrated Loop 10 µL		5/pk	5pk/cs	MP190-25
Calibrated Loop 1 µL		5/pk	5pk/cs	MP195-25
Small Loop	2.26mm	25/pk	—	MP197-25
Medium Loop	2.91mm	25/pk	—	MP199-25
Large Loop	3.57mm	25/pk	—	MP198-25

LeStreak®

Small Streaker Blue	2.26mm	5/pk	—	MP184-5
Medium Streaker Green	2.91mm	5/pk	—	MP182-5
Large Streaker Red	3.57mm	5/pk	—	MP180-5
Calibrated Streaker 1 µL		1 each	5/pk	MP185-5
Calibrated Streaker 10 µL		1 each	5/pk	MP187-5

LeStab® Needles

LeStab Micro Needle		5/pk	—	MP186-5
---------------------	--	------	---	---------



Platinum Inoculating Loops

Decon's Platinum Loops are high quality, durable loops designed for multiple re-use. The Platinum Loops are offered in four different loop dimensions. Handles sold separately. See below.

Platinum Inoculating Loops

Item	Loop ID	Shank Length	Std. Size	Alt. Size	Decon Cat #
Pt/15%Ir	2mm	45mm	1 each	10/Pk	MP 202-1
Pt/15%Ir	3mm	45mm	1 each	10/Pk	MP 203-1
Pt/3.5%Rh	4mm	25mm	1 each	10/Pk	MP 304-1
Pt/15%Ir	5mm	45mm	1 each	10/Pk	MP 205-1

LeLoop®/Platinum Loop Holders

Brass Loop Holder	6"		1 each	10/pk	MP6-10
Brass Loop Holder	8"		1 each	10/pk	MP8-10



HydeAway™

HydeAway aldehyde decontamination management products effectively convert used aldehydes into non-hazardous material to ensure safety and to facilitate disposal. HydeAway offers the highest neutralization ratio available and can be used for formaldehydes, formalin and gluteraldehydes. Crystalline Hyde-Away turns aldehydes into a solid polymer at a 1:30 ratio and the liquid version

turns aldehyde into a fluid suitable for sewer system disposal at a 1:15 ratio. Tested via EPA methods and proven safe for drain or trash disposal.

Item	Standard Size	Alternate Size	Decon Catalog #
Hyde-Away (Crystals)	32 oz.	12 x 32 oz.	2224
Hyde-Away (Liquid)	1 gallon (3.8L)	4 x 1 gallon	2201

Neutralize It!

HydeAway™ neutralizes formalin safely and completely...

- Neutralizes formaldehyde, gluteraldehyde
- Liquid for drain disposal
- Crystals for solidification
- Tested and safe for environment



Maximize It!

HydeAway™ has the best neutralization ratio by far...

- Liquid 1:15 - 1 gallon HydeAway neutralizes 15 gallons of formalin
- Crystals 1:30 - 32 oz HydeAway Neutralizes 7 gallons of formalin
- Closest competitors neutralization ratio: 1:5
- Save big \$ in neutralization disposal cost



Manufactured by Decon Laboratories, Inc.
Technical Information (800) 332-6647 or
www.deconlabs.com

Research Literature # BN1030039
Healthcare Literature # FHC-208
Canada Literature # DEC-040



460 Glennie Circle
King of Prussia, PA 19406
(800) 332-6647
www.deconlabs.com